


SOUND Veterinary Equipment

FEBRUARY - MARCH 2017

v-gel[®] ADVANCED VETERINARY AIRWAY MANAGEMENT SYSTEM

v-gel[®] is the first true veterinary species specific supraglottic airway device. Anatomically designed specifically for Rabbits and Cats. Manufactured from a soft gel-like material to give a high quality pressure seal that avoids laryngeal and tracheal trauma resulting in safer anaesthetics and superb comfortable patient recoveries.

- Fast, easy, safe and stress free device insertions.
- No post-operative coughing or gagging.
- Low airway breathing resistance due to the large airway channel.
- High quality pressure seal restricts leakage of agents.
- Super soft contoured tip for a highly effective upper oesophagus seal to prevent potential aspiration of reflux fluid.
- Integral bite block to stop patient damaging the device and occluding the airway.
- Designed to be reused and sterilised.

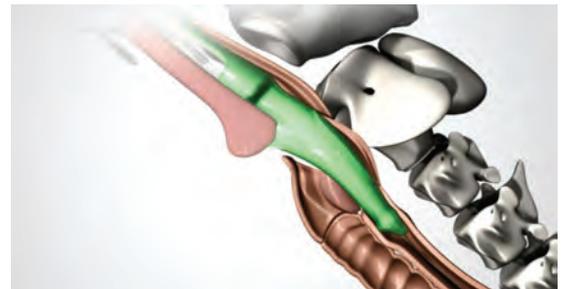

Rabbit v-gel Each	215⁰⁰
Cat v-gel Each	225⁰⁰
d-grip Tube Holder	145⁰⁰

Available in sizes 1-6, best selling sizes 2 & 3
Contact our office for sizing details



NEW

VETRONICS MINI TEMP

Hand-held battery operated temperature monitor. Can be used with standard YSI 400 series temperature probes. Maintain your patient in comfort zone 39.6°-36.6°C. Flashing LED indicates comfort zone. Accuracy 0.1°C. Veterinary specific.

MiniTemp 400 MV-250	290⁰⁰
MiniTemp Exotics MV-260	290⁰⁰



NEW

VETRONICS IMPACT III

Hand-held Capnograph with optional Pulse Oximetry and Temperature. Delivers performance and reliability for exotics and low flow applications. Micro-stream sampling for all animals over 100grams. Rechargeable battery operation for up to 6 hours. Contact us for further information.

Impact III Capnograph MV-110	3940⁰⁰
------------------------------	--------------------------

NEW



VETRONICS SAV04 VENTILATOR

Quality micro-processor controlled ventilator designed for small animals and exotics. Connects to fresh gas supply on anaesthesia machine. Measures and displays inspiratory time and respiratory rate. Battery back-up for emergency or mobile use. Clear & bright display.

SAV04 Ventilator MV-200	2260⁰⁰
-------------------------	--------------------------

WE NOW STOCK **Norfolk Vet Products** CONTACT OUR OFFICE

SOUND VETERINARY EQUIPMENT PTY LTD

UNIT 4, 19 CORNHILL ST, FERNTREE GULLY VIC 3156

P: 1300 881 681 | F: 1300 881 682 | E: contact@soundveterinary.com.au

Visit our E-Commerce Web Site: www.soundveterinary.com.au



GREEN CARPAL SPLINTS

Strong, lightweight plastic with foam padding. Splint may be contoured with moderate heat. Modified cutout for carpal pad. Green.

Small	2.5 x 20cm	3 ⁴⁰
Medium	3.5 x 30cm	3 ⁸⁵
Large	5 x 40cm	5 ⁰⁰



BLUE CARPAL SPLINTS

Made from a heavier material reinforced with longitudinal ribs. A recessed area reduces pressure on the carpal pad. A slight bend in the distal end follows the natural resting position. Open pore foam for comfort.

Toy	15 ⁹⁰	Large	20 ⁴⁰
Small	16 ⁸⁰	X-Large	26 ³⁰
Medium	17 ⁷⁰		



QUICK SPLINTS

Preformed Quick Splints made of high impact ABS plastic. Scored for individual sizing. Excellent replacement for Thomas Schroeder splints. Paired left and right sets.

FRONT		REAR	
Small	21 ⁸⁰	Small	18 ¹⁵
Medium	25 ⁸⁵	Medium	21 ⁵⁵
Large	31 ⁵⁰	Large	25 ⁸⁵



UV WOODS LAMP

Mains powered UV woods lamp with 5 diopter magnification lens (5 x 10cm). Four 4 watt UV light bulbs with bulb protection cover. Power cord length: 1.8 metres.

UV Woods Lamp
DL-900 **262⁴⁵**



CAGE CADDY

Cage caddy designed to hang on the cage door to store patient personal items such as collars and medication. Three separate compartments with lids. Blue polypropylene.

Size: 7.5 (H) x 32 (W) x 10 (D) cm

Cage Caddy
CA-100 **21⁷⁵**



RIGID STRETCHER

Rigid frame, 60 x 120cm. Lightweight and easy to clean. Cover is made of dirt resistant, durable nylon. Easily removed for machine washing. Secure with easy adjusting body straps with snap closures. Two wheels on one end.

Rigid Stretcher
AS-100 **131⁵⁰**



FERAL CAT OR SMALL ANIMAL DEN

Designed to fit into a standard cage this den provides a frightened or fractious animal a quiet hiding place. Reduces risk of injury to the handler while providing a stress free environment. Moulded plastic with carry handle, slide in front door and side opening porthole. 37 x 27 x 25cm.

Small Animal Den
AR-400 **213¹⁵**



CAT OR SMALL ANIMAL SHIELD

Just slip this light transparent panel into the cage, push toward the rear and inject through the evenly spaced holes. Easy to clean. Safe for the handler. Integral handle for one handed grip.

Size: 37.5 x 37.5cm.

Small Animal Shield
AR-300 **113⁴⁰**



CAT GRASPER

Lightweight grasper designed for effective, but safe control of cats. Ideal for reaching into the back of a cage to secure a fractious feline. Smooth covered jaws to minimise discomfort. If necessary the animal can be quickly released. 70cm long.

Cat Grasper
CG-100 **96⁰⁰**

ARTHREX VET SYSTEMS

TightRope CCL & FiberWire

TightRope CCL was developed to provide a minimally invasive and improved method for extracapsular stabilisation of the CCL. This technique does not require cutting of bone like the TPLO or TTA procedures. Instead, it uses small drill holes in the femur and tibia to pass a synthetic ligament-like biomaterial through a small incision to provide bone-to-bone stabilisation during healing. The biomaterial used for the TightRope CCL is called FiberTape. This is a Kevlarlike material that is used extensively in human surgery for many orthopaedic applications. This material has properties that make it stronger and less prone to failure than any other suture materials currently being used for CCL reconstructions. [Contact our office for additional information on the TightRope procedure or any additional products in the extensive Arthrex range.](#)



Standard TightRope Kit VAR-2800	335⁰⁰
Mini TightRope Kit VAR-2801	190⁰⁰
3.5mm Cannulated Drill Bit	160⁰⁰
For use with Standard TightRope VAR-8920DC	
2.7mm Cannulated Drill Bit	160⁰⁰
For use with Mini TightRope VAR-8911DC	
Guide Wires 0.049" Pack of 6 VAR-8920P	130⁰⁰



AUTOCLAVABLE ORTHOPAEDIC CORDLESS DRILL & SAW SETS

Fully autoclavable. Reliable German motors for smooth use and minimal noise. Cannulated drill takes up to 3mm pins while Jacobs Chuck capacity is 4.8mm. Range of blades available for Oscillating Saw.

Drill Set OI-790	3590⁰⁰
Saw Set OI-365	3590⁰⁰



BOSCH CORDLESS DRILL & SAW OSCILLATING SAW SETS

Bosch Professional cordless Drill & Saw. Supplied with 2 Lithium-ion batteries and autoclavable shrouds. Drill supplied with Jacobs Chuck and Drill Extension. Saw supplied with surgical converter.

Drill Set OI-780	617⁷⁰
Saw Set OI-360	745⁰⁰



SCREW RACKS WITH LID

Free standing screw racks with lid designed for a range of cortical & locking screws. Marked with sizes. 110 x 90 x 60mm (HxWxD)

1.5mm OIC-300	98⁰⁰
2.0mm OIC-301	98⁰⁰
2.4mm OIC-302	98⁰⁰
2.7mm OIC-303	98⁰⁰
3.5mm OIC-304	98⁰⁰



WELLIONVET GLUCOMETER

Veterinary specific blood glucose monitor – for reliability and accuracy. Easy operation with alarms. Large, clearly legible digits with back-lit display. Code chips specific for cats, dogs or horses.

Test Strips Box of 50 \$75⁰⁰ (BGM-210)

WellionVet Glucometer BGM-200	79⁵⁰
-------------------------------	------------------------



WELLIONVET BELUA BLOOD GLUCOSE & KETONES

Veterinary specific blood glucose monitor & ketone monitor specifically designed for cats, dogs and cows. Easy operation with alarms. Large, clearly legible digits with back-lit display.

Measure subclinical ketosis before visible symptoms appear.



BeluaVet Glucometer BGM-300	79⁵⁰
-----------------------------	------------------------

Blood Glucose Cats, Dogs & Cows Pk 50 BGM-310	75⁰⁰
Ketones Cats & Dogs Pk 10 BGM-320	60⁰⁰
Ketones Cows Pk 20 BGM-330	80⁰⁰



DISPOSABLE SURGICAL DRAPE MATERIAL

Disposable, blue, rich cloth-like feel of Sontara fabric packaged in a fan folded box. Simply cut to size and autoclave. Size 95cm x 91 metres.

Material 91mtr **118⁰⁰**
SDD-100



PORTABLE CAUTERY PEN

Portable Cautery Pen powered by 2 AA batteries. Full set includes four replaceable cautery tips. 1/2" Fine Electrode for hemostasis control, 1/2" Loop for wide use such as small skin tumours, 2" Fine Electrode for difficult areas & 2" Loop. High Temp.

Cautery Pen Set **106⁶⁰**
CU-200



BREATHE SAFE RESPIRATORY MONITOR

Great apnea alert monitor. Easy to use, just connects between e/t tube and breathing circuit. No adjustments needed, indicates every breath with a beep. A distinctive alert sounds if patient fails to breathe for 45 seconds. Auto shut-off. Long life lithium battery.

Breathe Safe **342⁸⁵**
MBS-100



MODIFIED JACKSON REES

Our most popular non-rebreathing system. The control valve allows you to close the system for hand bagging and then to reopen it to passive condition. Supplied with 0.5ltr bag.

Modified Jackson Rees **394⁵**
Non-rebreathing 32359B6



UNIVERSAL F CIRCUIT

Our most popular rebreathing system. Designed to maintain humidity and prevent heat loss. 1 metre in length, supplied with 2ltr breathing bag.

Universal F Circuit **394⁵**
Rebreathing 32359B3



K-9 EZ BOY POSITIONING AIDS

Suitable for radiographic or surgical positioning. V trough style, comfortable for the patient. Nylon reinforced vinyl. Features hook and loop straps to secure limbs. Easily cleaned for re-use.


Medium 37x17x12cm TP-210 **136⁰⁰**
Large 60x27x17cm TP-220 **179¹⁵**



D GRIP TUBE HOLDER

Fully adaptable airway management support system. Designed to support the load borne by anaesthetic circuits, airway devices, face masks, monitoring cables and sampling lines. Flexibility to adjust to the perfect position.


d-grip Tube Holder **145⁰⁰**
AVG-300

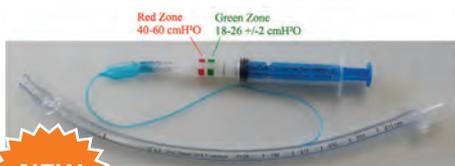

TRU-CUFF CUFF PRESSURE INDICATOR

The affordable inflation syringe with cuff pressure indicator. Correct cuff pressure is 20-30cm H₂O. Pressure above 30cm H₂O will exceed the surrounding tracheal blood supply pressure & constrict proper blood flow. Prolonged exposure to high pressure results in coughing and soreness. 3/Pack.

Tru-Cuff Syringe **35⁹⁰**
AET-TC (Pack of 3)

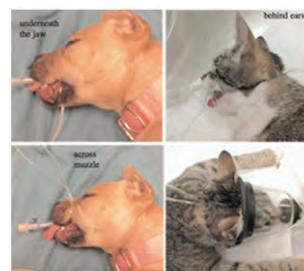

EZ TIES

Brilliant securing device! Ideal for endotracheal tubes, anaesthesia face masks, oxygen therapy masks etc. The EZ Tie adjusts for different size patients. Pack of 3.

EZ Ties **15⁴⁰**
AET-EZ

B.BRAUN VOLUMETRIC & SYRINGE INFUSION PUMPS

The B.Braun Perfusor Space Syringe Pump and Infusomat Space Series Infusion Pumps are setting new standards for veterinary infusion technology. Low running costs using standard giving set and unbeatable reliability. Made in Germany and backed by a 2 year warranty. We highly recommend these outstanding pumps.

- Identical user interface in both the infusion and syringe pump.
- Internal rechargeable battery.
- Supplied with pole mounting bracket.
- Pumps may be stacked on top of each other and locked together.
- Lighted keypad & display
- Networkable
- Higher accuracy compared with other pumps at low flow rates.
- Displays important information on the screen, visible at a glance.
- **Suitable for delivering IV Fluids and Blood.**
- Better error detection systems.
- Economical sets, from \$1.25
- 2 Year Warranty



IPBB-100	Infusomat Space P Infusion Pump	1675⁰⁰
IPBB-600	Intrafix Primline Set 20 Drops/ml	1²⁵
IPBB-700	Sangofix Blood Set 20 Drops/ml	2²⁵
IPBB-200	Perfusor Space Syringe Pump	2090⁰⁰

IM3 42-12 ULTRASONIC SCALER



Best settling 42-12 scaler. Adopts a gentle technique with rotational tip movement and ultrasonic removal of the plaque and calculus. Titanium tips.

42-12 Scaler
IM3-06979 **2600⁰⁰**

IM3 14 PIECE INSTRUMENT SET



General purpose dental set allowing you to examine, scale and extract teeth. Contains 14 instruments including Stubby Winged Elevators 1-5mm.

14 Piece Set
IM3-D1100 **395⁰⁰**

IM3 8 PIECE RABBIT & RODENT KIT



Specially designed Rabbit & Rodent Kit containing mouth gags, cheek dilators, extraction forceps, molar luxator, incisor luxator and spatula.

8 Piece Set
IM3-D1034 **295⁰⁰**

IM3 GS DELUXE & GS DELUXE LED DENTAL SYSTEMS

The iM3 GS Deluxe & iM3 GS Deluxe LED Dental Systems set the standard for dental machines. The unique iM3 design is your guarantee of many years of dependable service backed by the industry leading 5 Year Warranty on the control top and industry leading technical support.

Complete dental systems with advanced features not found on any other dental machine. Includes a quality push button high speed handpiece with water cooling for fast sectioning of teeth. The low speed handpiece and nose cone swivel for comfort and reduces strain on the hand. The low speed 4:1 handpiece is designed to use with iM3 R&R oscillating prophylaxis angles that won't wind up muzzle hair. A 3-way air water syringe with the ability to select either water or CLS. The GS Deluxe LED features the industry first self generating LED light. The fantastic performance of the LED will light up the mouth.



GS Deluxe
IM3-U3012 **6895⁰⁰**

GS Deluxe LED
IM3-U3019 **7395⁰⁰**

Order an iM3 Dental System with 42-12 Scaler and receive a FREE starter kit. Contact us for further information.

SURGIVET CDS9000 ANAESTHESIA MACHINE

5 Year Warranty

Outstanding quality, robust and reliable.

- Compact design
- Aluminium block with built-in inhalation & exhalation valves
- Maintenance free pop-off valve with pop-off occlusion valve
- Fresh gas port for non-rebreathing
- 200ml to 4ltr oxygen rotameter
- Easy, quick change absorber designed to use less soda-lime but allows expired gas to be exposed to more soda lime
- Suitable for patients up to 110kg
- Bench, wall or mobile models
- Vaporiser sold separately
- Supplied with all accessories including breathing circuits & bags, absorbent, scavenge & oxygen hoses.
- Outstanding quality machine with high end features.



Mobile CDS9000 V9000 **2780⁰⁰**

Bench Mount V9010 **2595⁰⁰**

Isoflurane Vaporiser Normally 1125⁰⁰ **995⁰⁰**

All CDS9000 Anaesthesia Machines are supplied with a bonus FREE monitor shelf valued at \$190. Contact our office for further information.

SUNTECH VET20 NIBP MONITOR

Quickly and easily provides reliable blood pressure measurements in a consult environment.



Features intuitive touch screen interface for viewing and analysing data. Measures systolic, diastolic, mean and HR. AA battery powered. 2 year warranty.

Suntech Vet20 MBP-300 **1850⁰⁰**

PETTRUST HAND HELD NIBP MONITOR

The PetTrust is a reliable and easy to use hand held blood pressure monitor. Ideal for spot checks or monitoring during surgery. Easy to use, displays systolic, diastolic and mean arterial pressure and heart rate. Cycle time feature.



PetTrust MBP-100 **1990⁰⁰**

SURGIVET ADVISOR VITAL SIGNS MONITOR

The leading choice for genuine veterinary specific monitoring

The SurgiVet Advisor features a large, easy to read display with intuitive design. The Advisor is available in multiple configurations and is also upgradable. Built in presets for the average cat, dog and horse. Effective and reliable monitoring for animals from 500grams to large animals.

A range of accessories are available to enhance the use of the monitor including a Mobile Roll Stand, Data Logger or Oesophageal ECG and Temperature Set.



Various configurations available. Contact our office for further information and pricing.

MODEL	ECG	HR	NIBP	SpO2	IBP	TEMP	ETCO2	RESP
V9230	✓	✓	✓	✓				
V9200	✓	✓	✓	✓				
V9200U	✓	✓	✓	✓		✓		
V9203	✓	✓	✓	✓	✓	✓	✓	✓